

E kangara aron anakin te **Rapid Antigen Test**

Angin mwaitin te 'Rapid Antigen Test' a iri aron tuangakia inano, ma tiaki ni kabane. Taiaoka ma ira nanon aia kaetieti taan karaoa te bai aei (manufacturer).

1
KATAURAOI

- Kakoaua ba ko teboka baim imwain ae ko waaki.
- Katauraoa te tabo ae itiaki ibukin karaoaan te tutuo – e raoiroi aon te taibora ae boraai ibukin karaoaan aei.
- Kawenei am bwai n tutuo – iai uoua te baeki (pouches).



2
SWAB

- Kauka te moan baeki (pouch) ao kaotinaoka
- the nasal swab (aekan taran te bai n ewa tanninga - ear bud)
- te ran ae uareke ae kawakinaki inanon ninnena (the sealed vial that has liquid in it), ao
- matan te vail ae uareke (the little vial cap).
- Karinnakoa te swab 2-3 cm inanon bangabangan te bairi are teuana – kaokoa nanon te bairi, ni katobibia nanon te bairi 4-5 te tai ao imwina kaotinaoka.
- Karinnakoa riki te swab anne inanon bangabangan te bairi are te kauoua ao katotonga are ko karaoia mai imwaina.



3
RIKOA AM TIAMBO (SAMPLE)

- Kauka te vial, karaurau n aki kabaroa te ran are inanova.
- Kateboa te swap inanon te ran anne (vial) nako nano ao kaminominoa inanon te ran anne ba e na mwaimwai ni kabane rabwatana.
- Kaotinaoka te swab, ma konna ni kanakoa te ran man te swab anne. E kakawaki ba e na kona n tinako ranin te swap nako nanon te vial.
- Kaina matan te vial ao ioioia ni kakibekibea man te iterana nakon te iterana ke ni kakibea n tabonibaim.
- Katikua te vial ni kateia ba ena tei raoi – ko kona ni kateia inanon te kirati (glass).



4
TUTUO

- Kauka te kauoua ni baeki (pouch) ao kaotinaoka te raurau n tutuo (testing tray).
- Ana am vial, kabarakia nako nano, ao katinakoa teutana 4 te timtim nako nanon te mronron (circular hole) iaon te raurau n tutuo.
- Tanninga inanon 15 te miniti.




5
MWIN REKENA / TE MWANEKA ARE IMWINA

Imwin 15 te miniti ao nora te tabo ae baoki i aon te raurau n tutuo

- Teuana te rain ae uraura iaan C = e aki kuneaki te manin aoraki (virus)
- Uoua te rain – uraura iaan te C ao te rain ae buru iaan T = e kuneaki te manin aoraki (virus) (e ngae ngkana e tane mauna te rain ae buru)
- Teuana te rain ae buru iaan C = e aki bainaki (invalid) mwin rekena

- Ribotina mwin rekena i aon Mycovidrecord.health.nz, ke tarebonia 0800 222 478 ao kotea are te ka-3 ni batin (option)
- Ngkana e kuneaki manin te aoraki n am tutuo:
 - taiaoka kaongoa am tia kam mwakuri ao koraki n aomata ae ko a tia n irengan (contact) ma ngaia.
 - ko na tangiria n tiku n onoti. Ibukin rongorongon te tiku n onoti ae boou nakon te itera ni bae n are Advice for people with COVID-19 i aon te Ministry of Health website.

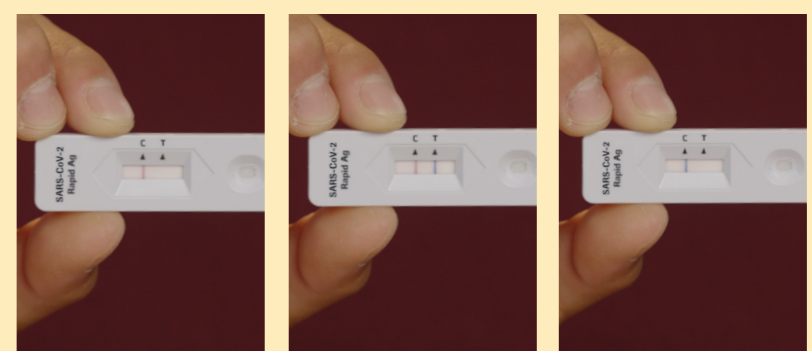

Not detected Detected Failed result