

# E fai pefea te faiga ote **Rapid Antigen Test**

Te lasiga o Rapid Antigen Test kits e masani o tautali ki fakatonuga mai lalo, kae sei ko mea katoa. Fakamolemole tautali ki fakatonuga o kamupane gaosi

## 1 FAKATOKA

- Masaua o fulu ou lima koi tuai koe o kaamata
- Fakatoka se koga lei mo fai tau test - iluga i te taipola mafola io me se teesi lei
- Fola fakalei a mea test a koe - e lua au taga



## 2 SWAB

- Tala te taga muamua o futi ki tua
  - te swab ote isu (e foliga pela mese mea fulu taliga)
  - te fagu e mau tena mata e isi ne vai iloto mote
  - pono ote fagu foliki
- Faulu te swab pela me 2-3 cm kiloto ite puisu e tasi - milimili iloto i tou isu, ko fakatamilo iei ite isu kaatoa 4-5 taimi ko tapale iei ki tua.
- Faulu te swab o fai foki pena i tou sua isu



## 3 TAE AU SWAB

- Tala te fagu fai fakalei moaa e moti se vai
- Faulu te swab kiloto ke tuu ki lalo ite muli ote fagu, oti ko fakatamilomilo iei ke miomio toe miomio ake mautinoa e pulutaki fakalei tau swab fakatasi mote vai
- Tapale te swab ki tua ko fakapiiki iei te mata i taimi e too mai ki tua. E taaua kii ke lasi te sua e kumi mai te swab ke fano kiloto ite fagu
- Fakapiki te pono iluga ite fagu ko lulu iei ki te sua feitu sua feitu io me petupetu ne koe ki ou maikao ke pulutaki fakalei
- Fakatafa pela te fagu kae tau loa o fakasaga te mata kiluga - e mafai o tuku ne koe iloto ite ipu kilasi




## 4 TEST

- Tala te lua o taga ko fakatafa keatea iei tau telee test
- Puke tau fagu, fuli te mata ki lalo, ko kumi malielie iei e 4 motiga ote vai kiloto ite puuga foliki
- Fakatali 15 minute



## 5 IKUGA/SUA STEP

- Mai tua ote 15 minute ko onoono koe kite koga tela e sikuea iluga ite telee test.
  - Tasi te laina kula mai lalo ote C = e seki iloa te vaelesi
  - Lua laina - laina kula mai lalo ote C ako te puluu mai lalo ote T = a koe ko maua ne te vaelesi (e ui te laina puluu e se matea lei)
  - Tasi te laina puluu mai lalo ote C = Te ikuga e se fakaoga
- Lipooti au ikuga iluga ite Mycovidrecord.health.nz io me vili 0800 222 478 ko paki iei Option 3
- Kafai tau test koe ko maua ne te vaelesi
  - Fakamolemole fakailoa ki tou pule mo niisi tino kola ne fesokotaki mo koe
  - Koe ko tau o fakataukatea. Fia iloa ne koe mea manakomia mafai ko fakataukatea koe fano kite *Advice for people with COVID-19* website ite laupepa iluga ite Ministry of Health.

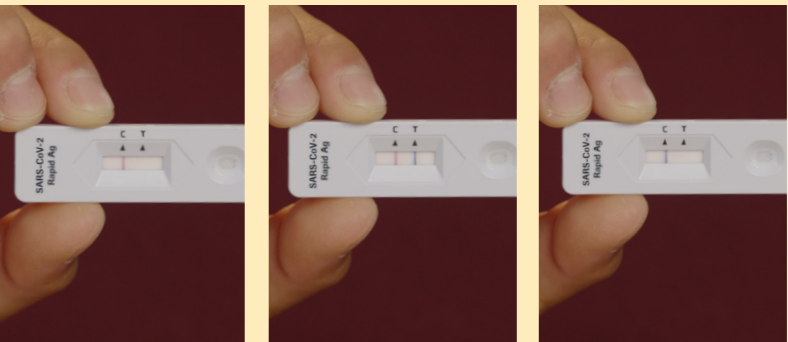

Not detected      Detected      Failed result